


1C53 Taleo Low Profile

Ready for everyday life.

Taleo Low Profile was designed for active individuals who navigate varied indoor and outdoor environments and place a high value on effortless walking and the ability to go wherever life takes them.

Designed for optimal outcomes:

- For users with limited clearance
- Dual springs and a long carbon base enable a smooth rollover and efficient energy return
- The flexible connection of the carbon springs in the forefoot, combined with multi-axial rotation movement result in optimal adaption to varying ground conditions
- The pyramid design allows a more natural ankle movement which provides a smooth rollover
- Three different heel wedges to optimize gait dynamics and reduce impact forces
- Protected against fresh, salt, and chlorine water
- Water runoff channels on connection adapter and drainage holes in footshell prevent water from collecting in the interior of the prosthesis

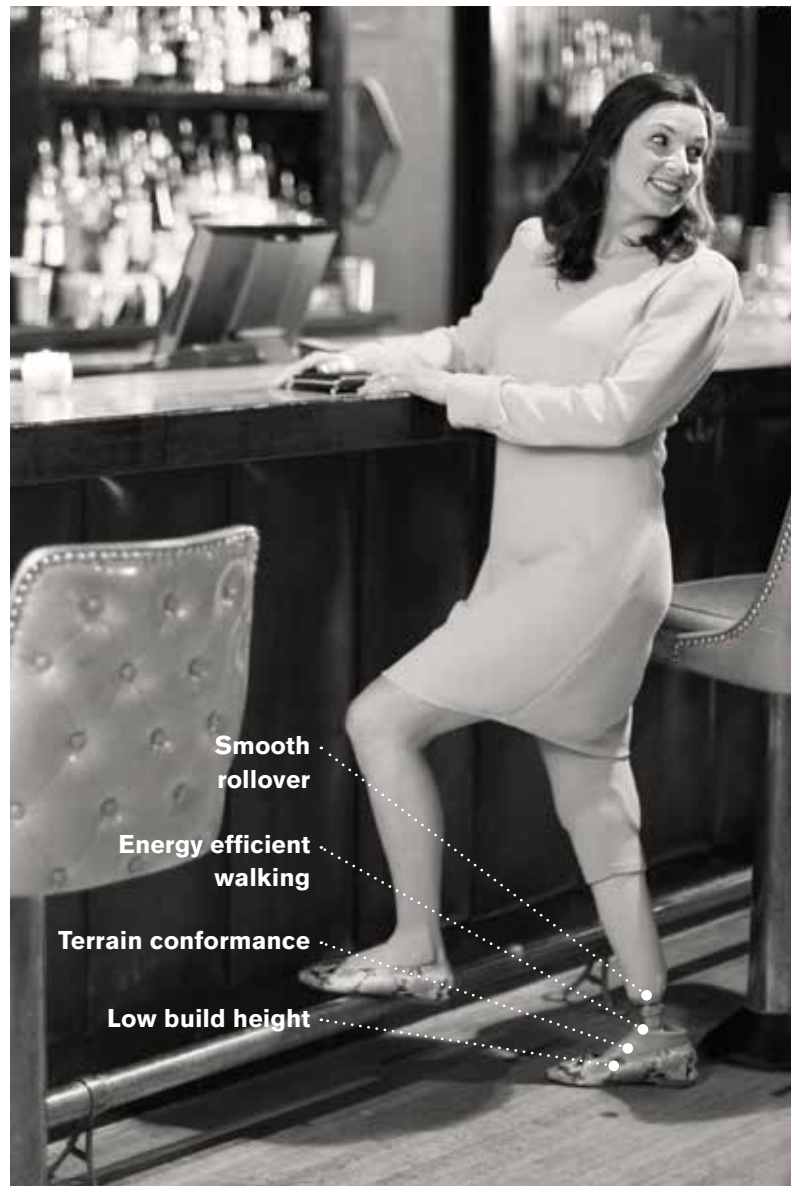

More than a foot. A foundation.

Technical Data

K Level	K3
Max. body weight	330 lbs / 150 kg
Sizes	22-30 cm
Footshell	Slim shape with 15 ± 5 mm heel height (size 22-25 cm) Normal shape with 10 ± 5 mm heel height (size 24-30 cm)
Weight without footshell*	12.52 oz / 355 g
Weight with normal footshell*	20.6 oz / 584 g
Build height	2 1/2 in / 65 mm (with normal footshell*)

*Reference size is 26 cm.



1C53 Taleo Low Profile



Order form


Customer		Shipping address (if different from customer address)	
Customer no.	<input type="text"/>	Customer no.	<input type="text"/>
Company	<input type="text"/>	Company	<input type="text"/>
Street	<input type="text"/>	Street	<input type="text"/>
Postal code/city	<input type="text"/>	Postal code/city	<input type="text"/>
O&P professional	<input type="text"/>	PO no.	<input type="text"/>

Selection of the spring stiffness relative to

1 body weight and activity as well as

Body weight lbs [kg]	Normal activity level	High activity level
Up to 113 [51]	1	2
114 – 129 [52 – 58]	2	3
130 – 149 [59 – 67]	3	4
150 – 170 [68 – 77]	4	5
171 – 195 [78 – 88]	5	6
196 – 220 [89 – 100]	6	7
221 – 254 [101 – 115]	7	8
255 – 287 [116 – 130]	8	9
288 – 330 [131 – 150]	9	–

2 foot size

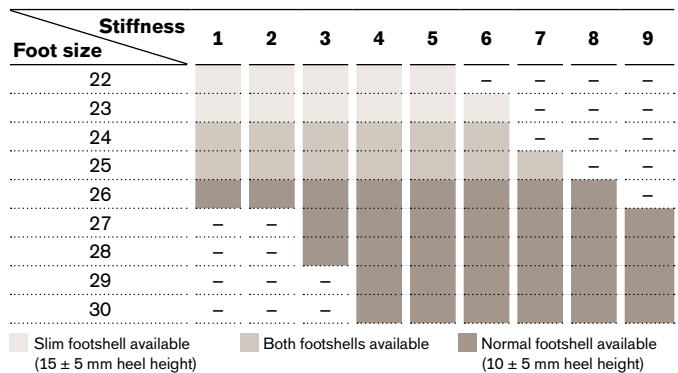

Example for ordering 1C53=R26-4-P/4N

Article no.	Side	Size	Stiffness	P /	Color	Shape
1C53=				P /		

Side	Size [cm]	Color	Shape
Right R	22, 23, ..., 30	Beige 4	Slim 22–25 cm
Left L		Light brown 15	Normal 24–30 cm



The scope of delivery includes the 1C53 foot module, the 2C15 footshell including connection cover, the 2F50 heel wedge set (comprising three different degrees of hardness) and a black Spectra sock.

Suggested combinations**

TF	KD	TT
7T452=2 SiOCX TF Seal & 6Y81 ProSeal liner Laminated socket & 6Y110 Skeo Sealing & 21Y21 ClickValve	Laminated socket & 6Y110 Skeo Sealing & 21Y21 ClickValve	4R180 Harmony P4 & 6Y512 Uneo 3D & 453A30/40 ProFlex Plus knee sleeve
3B1* Genium 3C98-3* C-Leg 4	3B1*=ST Genium 3C88-3* C-Leg 4	Laminated socket & 6Y92 or Unique custom liner & 453A30/40 ProFlex Plus knee sleeve & 4R136=EL one-way valve
4R57/4R57=ST Rotation adapter 4R57=WR/4R57=WR-ST Waterproof rotation adapter 4R10=110 Quickchange		4R10=111 Quickchange
3F1=2 + 99B120=* Functional Cosmesis for Genium 3F1=1 + 99B120=* Functional Cosmesis for C-Leg 4	3F1=2 + 99B120=* Functional Cosmesis for Genium 3F1=1 + 99B120=* Functional Cosmesis for C-Leg 4	

** All components are sold separately and are available Ottobock products that are compatible with the 1C53 Taleo Low Profile, which helps ensure optimal performance. Practitioners need to select components based upon individual patient criteria.